



Philips Healthcare

- 1/3 -

FSN86100207A

October 2019

**URGENT – Medical Device Correction**  
**Efficia DFM100 (Model number 866199)**

<Hospital/Clinic Name>

cc: Chairman Medical Board and relevant Head of Departments

Dear Customer,

Philips has identified that the Efficia DFM100 Defibrillator/Monitor (Model number 866199) rotary therapy selector switch may fail, resulting in unpredictable device behaviour. These behaviours include:

- The device may not turn on
- The device may not perform the selected function
- The device may deliver a shock with an energy level different than the setting selected by the user

Should one of these behaviours occur, appropriate therapy delivery may be delayed. To date, Philips has not received any reports of deaths resulting from this switch failure.

As a remedy, Philips will install a replacement switch in affected devices at no charge to the customer.

The purpose of this notification is to:

- Describe actions that you should take to mitigate risk to patients
- Recommend that unit be removed from service if they exhibit these symptoms
- Describe the corrective action planned by Philips to address the issue

**This document contains important information for the continued safe and proper use of your equipment**


Please review the following information with all members of your staff who should be aware of the contents of this communication.

Please retain a copy with the equipment Instructions for Use.

Follow the "ACTION TO BE TAKEN BY CUSTOMER / USER" section of the notice.

If you have questions regarding this notification or need any further information or support, please contact your local Philips representative. <Philips representative contact details to be completed by the KM / country>, Healthcare Representative/Modality Engineer: 1800-744-5477 or (Overseas Number).

Sincerely,

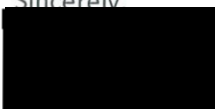


Gregory M Ayers, MD, PhD  
Head of Post Market Surveillance  
Associate Chief Medical Officer  
Monitoring & Analytics and Therapeutic Care



Li Ping  
Senior Q&R Manager, MA&TC Q&R PQMS

## URGENT – Medical Device Correction Efficia DFM100 (Model number 866199)
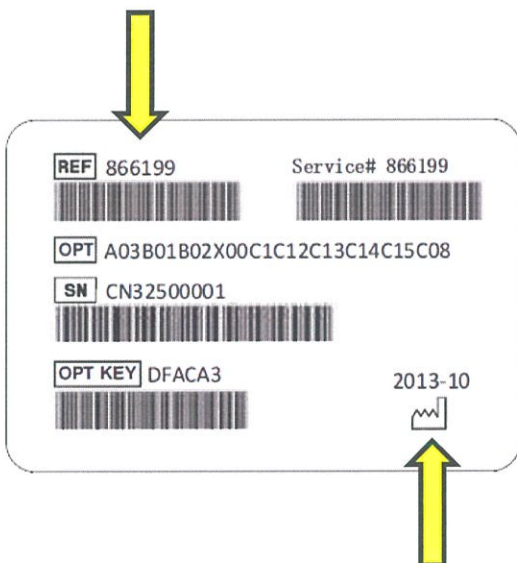
<b>AFFECTED PRODUCTS</b>	All Philips Efficia DFM100 Defibrillator/Monitor manufactured (Model number 866199) prior to 1 May 2017
<b>HOW TO IDENTIFY AFFECTED PRODUCTS</b>	<p>The model number of the Philips Efficia DFM100 is printed on the primary label on the bottom of the device</p>  <p>The date of manufacture of the Philips Efficia DFM100 is printed on the primary label on the bottom of the device</p>
<b>BEHAVIOUR DESCRIPTION</b>	<p>The Philips Efficia DFM100 Defibrillator/Monitor rotary therapy selector switch may fail, resulting in unexpected device behaviour. These behaviours include:</p> <ul style="list-style-type: none"> <li>• The device may not turn on</li> <li>• The device may not perform the selected function</li> <li>• The device may deliver a shock with an energy level different from the setting selected by the user</li> </ul>
<b>HAZARD INVOLVED</b>	<p>These device behaviours could result in a delay in therapy or failure to deliver the intended therapy.</p> <p>Philips has not received any reports of patient deaths associated with this failure of an Efficia DFM100 Monitor/Defibrillator.</p>



**URGENT – Medical Device Correction**  
**Efficia DFM100 (Model number 866199)**

<b>ACTION TO BE TAKEN BY CUSTOMER / USER</b>	<p>The device is safe to use and can remain in service if it does not exhibit any of these behaviors described in this Notice.</p> <p>Continue to perform Shift Checks and Operational checks as recommended in the Instructions for Use (IFU) as this reduces the risk of a failure during use.</p> <p>If you identify a device that exhibit any of these behaviors, please remove it from service and contact Philips to request service.</p>
<b>ACTIONS PLANNED BY PHILIPS</b>	<p>Philips will contact you to arrange for repair of your unit. Philips will install a replacement switch in affected devices at no charge to the customer.</p>
<b>FURTHER INFORMATION AND SUPPORT</b>	<p>If you need further information or support concerning this notification, please contact your local Philips representative &lt;Healthcare Representative /Modality Engineer: 1800-744-5477 or (Overseas Number)&gt;.</p>